

Product Manual

Revision date: 15 Oct 2019

Serotonin / 5-Hydroxytryptamine (5-HT) ELISA Kit

Catalog No.: BTA14693

Size: 96T

Range: 15.63 ng/ml - 1000 ng/ml

Sensitivity: 9.38 ng/ml

Storage: Store the 96-well plate, Standards, HRP-conjugate reagent and Biotin-conjugated antibody at -20°C, and the rest of the kit components at 4°C.

Application: For quantitative detection of 5-HT in Serum, Plasma and other biological fluids.

Introduction: Serotonin or 5-hydroxytryptamine (5-HT) is a monoamine neurotransmitter. Biochemically derived from tryptophan, serotonin is primarily found in the gastrointestinal tract (GI tract), blood platelets, and the central nervous system (CNS) of animals, including humans. It is popularly thought to be a contributor to feelings of well-being and happiness. Serotonin is a neurotransmitter and is found in all bilateral animals, where it mediates gut movements and the animal's perceptions of resource availability. In less complex animals, such as some invertebrates, resources simply mean food availability. In more complex animals, such as arthropods and vertebrates, resources also can mean social dominance. In response to the perceived abundance or scarcity of resources, an animal's growth, reproduction or mood may be elevated or lowered. This may somewhat depend on how much serotonin the organism has at its disposal.

Principle of the Assay

This kit is based on a competitive enzyme-linked immuno-sorbent assay technology. 5-HT is pre-coated onto a 96-well plate. The standards, samples and a biotin conjugated antibody specific to 5-HT are added to the wells and incubated. After washing away the unbound conjugates, avidin conjugated to Horseradish Peroxidase is added to each microplate well. After TMB substrate solution is added only wells that contain 5-HT will produce a blue colour product that changes into yellow after adding the stop solution. The intensity of the yellow colour is inversely proportional to the 5-HT amount bound on the plate. The O.D. absorbance is measured spectrophotometrically at 450 nm in a microplate reader, and then the concentration of 5-HT can be calculated.

Kit components

- | | |
|---|--|
| 1. One pre-coated 96-well microplate (12 × 8 well strips) | 4. tips |
| 2. Standard: 2 vials | 5. Squirt bottle or automated microplate washer |
| 3. Sample/Standard Diluent Buffer: 20 ml | 6. ELISA shaker |
| 4. Biotin conjugated antibody (100X): 120 µl | 7. 1.5 ml tubes to prepare standard/sample dilutions |
| 5. Antibody diluent buffer: 12 ml | 8. Distilled or deionized water |
| 6. HRP-conjugate reagent (100X): 120 µl | 9. Absorbent filter papers |
| 7. HRP diluent buffer: 12 ml | 10. 100 ml and 1 liter graduated cylinders |
| 8. TMB substrate reagent: 10 ml | |
| 9. Stop solution: 10 ml | |
| 10. Wash buffer (25X): 30 ml | |
| 11. Plate sealer: 5 pieces | |

Material Required But Not Provided

- 37°C incubator
- Microplate reader (wavelength: 450 nm)
- Multi and single channel pipettes and sterile pipette

Product Manual

Revision date: 15 Oct 2019

Protocol

A. Preparation of sample and reagents

A. Sample

Store samples to be assayed within 24 hours at 2-8°C. Alternatively, aliquot and store at -20°C or -80°C for long term. Avoid

- **Serum:** should be collected into a serum separator tube. Coagulate the serum by leaving the tube undisturbed in a vertical position overnight at 4°C or at room temperature for up to 60 minutes. Centrifuge at approximately 1000 × g for 20 min. Analyze the serum immediately or aliquot and store at -20°C or -80°C.
- **Plasma:** plasma using heparin or EDTA as an anticoagulant. Centrifuge for 15 minutes at 1000 × g within 30 minutes of collection. Assay immediately or aliquot and store at -20°C. Avoid hemolysis and high cholesterol samples.
- **Other biological fluids:** at approximately 1000 × g for 20 min to remove precipitant. Analyze immediately or aliquot and store at -20°C or -80°C.

Note:

- » Fresh samples or recently obtained samples are recommended to prevent degradation and denaturalization that may lead to erroneous results. It is recommended to store samples to be used within 5 days at 4°C, within 1 month at -20°C and within 2 months at -80°C.
- » Samples should be clear and transparent. Samples must be diluted so that the expected concentration falls within the kit's range.
- » Please bring sample slowly to room temperature. Sample hemolysis will influence the result. Hemolyzed specimen should not be used. Samples that contain NaN₃ cannot be detected as it interferes with HRP.
- » Always use non-pyrogenic, endotoxin-free tubes for blood collection.


1. Wash buffer

Dilute the concentrated Wash buffer 25-fold (1/25) with distilled water (i.e. add 15 ml of concentrated wash buffer into 360 ml of distilled water).

1. Standard

Preparation of the 5-HT standard: standard solution should be prepared no more than 15 min prior to the experiment. Centrifuge at 10,000×g for 1 minute.

2. 1000 ng/ml standard solution. Add 1 ml of Sample/Standard diluent buffer into one Standard tube. Allow the reconstituted standard to sit for 15 minutes with gentle agitation prior to carrying out the serial dilutions; avoiding foaming or bubbles.
3. 500 ng/ml → 15.625 ng/ml standard solutions: Label 6 tubes with 500 ng/ml, 250 ng/ml, 125 ng/ml, 62.5 ng/ml, 31.25 ng/ml and 15.625 ng/ml. Aliquot 0.5 ml of the Sample / Standard diluent buffer into each tube. Add 0.5 ml of the above 1000 ng/ml standard solution into 1st tube and mix thoroughly. Transfer 0.5 ml from 1st tube to 2nd tube and mix thoroughly. Transfer 0.5 ml from 2nd tube to 3rd tube and mix thoroughly, and so on.

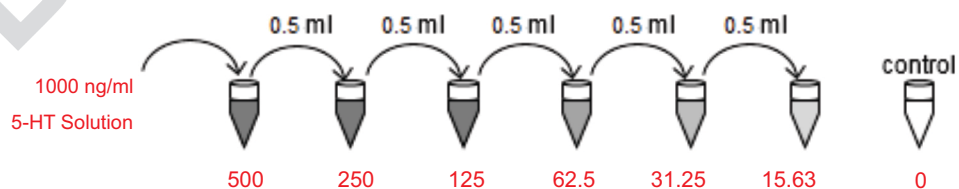

Note: Do not vortex the standard during reconstitution, as this will destabilize the protein. Once your standard has been reconstituted, it should be used right away. Use the diluted Standards for a single assay procedure and discard after use.

1. **Preparation of Biotin Detection Reagent working solution:** prepare no more than 1 hour before the experiment.
2. Calculate the total volume of the working solution: 0.05 ml / well × quantity of wells. (Allow 0.1-0.2 ml more than the total volume).

Product Manual

Revision date: 15 Oct 2019

- a.) Dilute the Biotin Detection Reagent with Detection Reagent diluent buffer at 1/100 and mix thoroughly. i.e. Add 1 μ l of Biotin Detection Reagent into 99 μ l of Detection Reagent diluent buffer. Discard after use.
- b.) **Preparation of HRP Conjugated Reagent working solution:** prepare no more than 30 min. before the experiment
- c.) Calculate the total volume of the working solution: 0.1 ml / well \times quantity of wells. (Allow 0.1-0.2 ml more than the total volume).
- d.) Dilute the HRP Conjugate Reagent with HRP diluent buffer at 1/100 and mix thoroughly. i.e. Add 1 μ l of HRP Conjugate Reagent into 99 μ l of HRP diluent buffer. Discard after use.

- a.) **Assay Procedure**

Equilibrate the kit components and samples to room temperature before use. It is recommended to plot a standard curve for each test.

- b.) **Set standard, test sample and control (zero) wells on the pre-coated plate and record their positions. It is recommended to measure each standard and sample in duplicate. Add the solution at the bottom of each well without touching the side walls. Mix the standards and samples up and down to be homogeneous before adding into the wells but avoid adding bubbles.**
- c.) Add 50 μ l of the prepared standards solutions into the standard wells. Add 50 μ l of Sample / Standard diluent buffer into the control (zero) wells.
- d.) Add 50 μ l of appropriately diluted sample into test sample wells.
- e.) Immediately add 50 μ l of Biotin conjugated antibody working solution into each well. Add the solution at the bottom of each well without touching the side wall.
- f.) Cover the plate with the plate sealer. Gently tap the plate to mix thoroughly. Incubate at 37°C for 45 minutes.
- g.) Remove the cover and discard the solution. Wash the plate 3 times with 1X Wash Buffer. Fill each well completely with Wash Buffer (350 μ l) using a multi-channel pipette or autowasher (1-2 minute soaking period is recommended). Complete removal of liquid at each step is essential for good performance. After the final wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean absorbent paper towels.
- h.) Add 100 μ l of HRP-conjugate working solution into each well. Cover the plate with a new sealer and incubate at 37°C for 30 minutes.
- i.) Remove the cover, discard the liquid and wash the plate 5 times with Wash buffer as explained in step 6.
- j.) Add 90 μ l of TMB Substrate into each well. Seal the plate with a cover and incubate at 37°C for 10-20 min. Avoid exposure to light. The incubation time is for reference only, the optimal time should be determined by end user. Do not exceed 30 min.
- k.) Add 50 μ l of Stop solution into each well to stop the enzyme reaction. It is important that the Stop Solution is mixed quickly and uniformly throughout the microplate to inactivate the enzyme completely.
- l.) Ensure that there are no fingerprints or water on the bottom of the plate, and that the fluid in the wells is free of bubbles. Measure the absorbance at 450 nm immediately.

This assay is competitive, therefore there is an inverse correlation between 5-HT concentration in the sample and the absorbance measured. Create a graph with the log of the standard concentration (y-axis) and absorbance measured (x-axis). Apply a best fit trendline through the standard points. Use this graph calculate sample concentrations based on their OD values. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

Note: If the samples measured were diluted, multiply the dilution factor by the concentrations from interpolation to obtain the concentration before dilution.

Product Manual

Revision date: 15 Oct 2019

A. Precautions

- B. Before using the kit, centrifuge the tubes briefly to bring down the contents trapped in the lid.
- C. If crystals have formed in the concentrated Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.
- D. Avoid foaming or bubbles when mixing or reconstituting components. Prepare the Standard dilutions within 15 min of use and discard any unused working standards. For each step in the procedure, total dispensing time for addition of reagents to the assay plate should not exceed 10 minutes.
- E. Do not let the wells uncovered for extended periods between incubation. Once reagents are added to the wells, avoid letting the strips dry as this can inactivate the biological material on the plate. Incubation time and temperature must be controlled.
- F. Ensure plates are properly sealed or covered during incubation steps.
- G. Complete removal of all solutions and buffers during wash steps is necessary for accurate measurement readings.
- H. Do not reuse pipette tips and tubes to avoid cross contamination.
- I. The TMB Substrate solution is easily contaminated; work under sterile conditions when handling the TMB substrate solution. The TMB Substrate solution should also be protected from light. Unreacted substrate should be colorless or very light yellow in appearance. Aspirate the dosage needed with sterilized tips and do not dump the residual solution back into the vial.

J. Precision

Intra-assay Precision (Precision within an assay): 3 samples with low, medium and high levels of 5-HT were tested 20 times on one plate, respectively.

Inter-assay Precision (Precision between assays): 3 samples with low, medium and high levels of 5-HT were tested on 3 different plates, 8 replicates in each plate.

CV (%) = (Standard Deviation / mean) × 100

Intra-Assay CV: ≤ 10%

Inter-Assay CV: ≤ 10%

Product Manual

Revision date: 15 Oct 2019

A. Typical Data & Standard Curve

Typical Standard Curve Data provided for demonstration purposes only. A new standard curve must be generated for each assay performed.

Concentration ng/ml OD450	0	15.625	31.25	62.5	125	250	500	1000
	2.451	1.735	1.409	1.120	0.834	0.555	0.372	0.266

